

SN: 08/095,555
FD: 7-26-93

S-73,893

PATENTS-US--A8095555

W-31-109-ENG-38

Philip E. Horwitz
Mark L. Dietz

PROCESS FOR THE SEPARATION
AND PURIFICATION OF
YTTRIUM-90 FOR MEDICAL
APPLICATIONS

095,555

DISCLAIMER

This report was prepared as an account of work sponsored by an agency of the United States Government. Neither the United States Government nor any agency thereof, nor any of their employees, makes any warranty, express or implied, or assumes any legal liability or responsibility for the accuracy, completeness, or usefulness of any information, apparatus, product, or process disclosed, or represents that its use would not infringe privately owned rights. Reference herein to any specific commercial product, process, or service by trade name, trademark, manufacturer, or otherwise does not necessarily constitute or imply its endorsement, recommendation, or favoring by the United States Government or any agency thereof. The views and opinions of authors expressed herein do not necessarily state or reflect those of the United States Government or any agency thereof.

A PROCESS FOR THE SEPARATION AND PURIFICATION OF YTTRIUM-90 FOR MEDICAL APPLICATIONS

Philip E. Horwitz
Mark L. Dietz

MASTER

DISTRIBUTION RESTRICTED TO U.S. ONLY

10

DISCLAIMER

Portions of this document may be illegible in electronic image products. Images are produced from the best available original document.

CONTRACTUAL ORIGIN OF THE INVENTION

The United States Government has rights in this invention pursuant to Contract No. W-31-109-ENG-38 between the U.S. Department of Energy and the University of Chicago.

5

BACKGROUND OF THE INVENTION


This invention relates to a process for the separation and purification of curie quantities of yttrium-90, produced from strontium-90 for use in the treatment of cancer and rheumatoid arthritis of the knee joint. More specifically, the invention relates to a process for the separation and purification of yttrium-90 while preserving the purity of the strontium-90 source, and removing other impurities from the yttrium-90.

10

5 Considerable effort has recently been directed toward the development
of site-specific methods for the treatment of various forms of cancer using
radionuclides. The suitability of a given radionuclide for use in such an
application is determined by several factors, among them its mode of
10 production, the availability of effective methods for its attachment to a
site-specific agent, its potential toxicity if detached from the agent, and its
therapeutic effectiveness. These considerations, together with the requirement
that the isotope have an appropriate half-life and decay scheme, severely limit
the number of radionuclides which may be seriously considered for use in
15 radioimmunotherapy.

Among the more attractive radionuclides for therapeutic applications is
yttrium-90. Its relatively short half-life (64.0 h) and maximum beta energy
(2.28 MeV) make it well suited for a variety of applications, ranging from the
radiolabeling of antibodies for tumor therapy to the production of
15 radiolabeled particles for the treatment of liver malignancies.

Yttrium-90 results from the decay of strontium-90 according to the following
scheme:

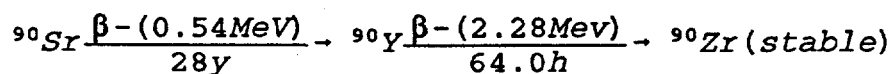

Before it can be safely employed in clinical applications, the ${}^{90}\text{Y}$ must
be made essentially free of ${}^{90}\text{Sr}$, an isotope known to cause bone marrow
20 suppression. In addition, any trace elements which could interfere with the

radiolabeling process by competing with ^{90}Y for binding sites must be removed.

5 Many methods for effecting the separation of ^{90}Y and ^{90}Sr have been described, among them solvent extraction, ion-exchange, precipitation, and various forms of chromatography. Of these, ion-exchange methods have probably received the most attention. Numerous procedures have been reported, for example, in which a cation exchange resin (e.g. Dowex 50) is used to retain ^{90}Sr , while the ^{90}Y is eluted with an aqueous complexant such as lactate, acetate, citrate, oxalate, or EDTA. Several of these procedures have been proposed as the basis for ^{90}Y generator systems. Unfortunately, in each case, the ^{90}Y is not eluted in a form suitable for direct labeling of antibodies. That is, often the concentration of the complexing agent is such that it will compete effectively with binding sites on the antibody for the activity, resulting in a decrease in the labeling efficiency. Thus, it usually becomes necessary to either remove these materials prior to antibody labeling or to carry out postlabeling purification. In addition, ion-exchange resins are plagued by a gradual loss of capacity due to radiation damage. As a result, ion-exchange is considered suitable only for subcurie quantities of activity, less than the quantities which often must be processed for clinical applications. Finally, achieving acceptable ^{90}Y yields while minimizing ^{90}Sr breakthrough often requires the use of long ion-exchange columns and impractically large volumes of eluent.

Kanapilly and Newton (1971) have described a process applicable to the separation of even multi-Ci quantities of ^{90}Y from ^{90}Sr which overcomes many of the limitations of ion-exchange methods. This process, which involves the precipitation of ^{90}Y as the phosphate, however, requires the addition of nonradioactive yttrium as a carrier, yielding a ^{90}Y product which is not, obviously, carrier free. More recently, workers at Oak Ridge National Laboratory have introduced a procedure involving an initial solvent extraction of the ^{90}Y from a dilute acid solution of $^{90}\text{Sr}/^{90}\text{Y}$ using bis (2-ethylhexyl) phosphoric acid in dodecane. Although good decontamination factors are routinely obtained, the useful life of the generator is limited, as radiolysis products of the extractant employed gradually accumulate in the ^{90}Sr stock. Moreover, the complexity of the process, which involves repeated stripping of the initial extractant solution to reduce trace impurities and repeated wet ashing of the stock solution to destroy dissolved organic phosphates, is prohibitive.

U.S. Patent 5,100,585, and U.S. Application Serial No. _____ (S-72,219) filed June 15, 1993, both assigned to the U.S. Department of Energy, and incorporated herein by reference, describe processes for the recovery of strontium and technetium values from acidic feed solutions containing other fission product values and containing up to 6 molar nitric acid. Recent work patented by the inventors, directed at improving analytical methodology for the determination of actinides and

radiostrontium in biological and environmental samples, has led to the development of extraction chromatographic resins capable of selectively sorbing actinides and lanthanides or radiostrontium from nitric acid solution. See, Horwitz et. al, "Method for Liquid Chromatographic Extraction of Strontium from Acid Solutions", U.S. Pat. No. 5,110,474, and Horwitz et. al, "Method for the Concentration and Separation of Actinides from Biological and Environmental Samples", U.S. Pat. No. 4,835,107. The present application discloses the application of these materials to the preparation of ^{90}Y of sufficient chemical and radiochemical purity to be suitable for use in medical applications.

Accordingly, it is an object of the present invention to provide a process for the preparation of ^{90}Y of sufficient chemical and radiochemical purity to be suitable for use in medical applications.

It is another object of the present invention to provide a process for the separation of ^{90}Y from ^{90}Sr by selectively separating the ^{90}Sr from the ^{90}Y .

Yet another object of the present invention is to provide a process for the separation of ^{90}Y from ^{90}Sr which avoids the need to frequently replace the ^{90}Sr and the generation of large quantities of radioactive waste.

SUMMARY OF THE INVENTION

The process for the separation and purification of yttrium-90 from a stock solution of strontium-90 involves an initial purification of the ^{90}Sr stock solution by dissolving the residue remaining after its evaporation in nitric acid

and passing it through an extraction chromatographic column containing a strontium-selective extractant sorbed on an inert polymeric support.

Strontium is retained, while decay products and impurities pass through and are discarded. The strontium is then removed from the column with dilute
5 nitric acid. Any traces of organic materials present in this purified ^{90}Sr are then removed by passing the solution through a bed of a polymeric sorbent.

The ^{90}Sr solution is then stored for a sufficient period of time to allow ^{90}Y ingrowth to occur in the solution, thus providing a $^{90}\text{Sr}/^{90}\text{Y}$ mixture. Next, the $^{90}\text{Sr}/^{90}\text{Y}$ mixture is acidified and passed through a series of strontium-

10 selective chromatographic columns such that essentially all of the ^{90}Sr is retained on the columns and the effluent from the third column contains essentially only ^{90}Y . The ^{90}Y effluent is then passed through a rare-earth selective chromatographic column such that ^{90}Y is retained, while any residual ^{90}Sr is stripped away. Sorbed ^{90}Y is then eluted from the rare-earth selective
15 column. Finally, each of the strontium-selective extractant columns is stripped of sorbed ^{90}Sr with dilute nitric acid and returning to storage to allow for new ^{90}Y growth. The strontium selective extractant resin may be obtained from an octanol solution of di-*t*-butylcyclohexano-18-crown-6, while the yttrium-selective extractant may be obtained from a solution of octyl(phenyl)-*N,N*-
20 diisobutylcarbamoylmethylphosphine oxide in tributylphosphate (CMPO in TBP). The various stages of the process are each physically isolated from each other and carried out in individual containment boxes. The total

strontium loaded onto the first column typically does not exceed twenty percent of its capacity.

BRIEF DESCRIPTION OF THE DRAWINGS

5 The above-mentioned and other features of the invention will become more apparent and be best understood, together with the description, by reference to the accompanying drawings, in which:

Figure 1 shows a column arrangement for the isolation of yttrium-90 from strontium-90 stock in accordance with the present invention;

10 Figure 2 shows a graphical representation of nitric acid dependencies of the retention of yttrium (III) and iron (III) on the rare-earth Selective chromatographic resin;

Figure 3 shows the effect of resin irradiation on the elution behavior of yttrium-90 on a rare-earth selective column;

15 Figure 4 shows the effect of resin irradiation on the stripping of strontium from a strontium-selective column; and,

Figure 5 shows the effect of resin irradiation on the retention of strontium by a strontium-selective column.

DESCRIPTION OF THE PREFERRED EMBODIMENTS

20 Figure 1 depicts the basic column arrangement and sequence of steps required to isolate ^{90}Y from the original ^{90}Sr stock solution. First, a preliminary purification of the ^{90}Sr stock is carried out by dissolving the residue remaining after evaporation of the original solution in 3 M nitric acid

and passing it through a chromatographic a column containing a strontium-selective extraction chromatographic resin. This resin is prepared with purified Amberchrom CG-71ms resin impregnated with an octanol solution of di-t-butylcyclohexano-18-crown-6 (hereafter abbreviated as DtBuCH18C6).

5 The resultant material, referred to hereafter as Sr•SpecTM column material (for strontium specific), is commercially available from ElChroM Industries, Inc. (Darien, IL). Strontium is retained, while any yttrium and zirconium which have accumulated via decay, as well as any impurities initially present, are flushed from the column with several free column volumes of 3 M nitric
10 acid (as shown in Fig. 1) and discarded. The purified strontium is then stripped from the column with dilute (~ 0.05 M) nitric acid and the solution passed through a bed of uncoated Amberlite XAD-7 resin (Rohm & Haas, Philadelphia, PA) to remove any traces of organic materials present. The treated solution is then either stored without further treatment or evaporated
15 to dryness. After a suitable period of ⁹⁰Y ingrowth, the ⁹⁰Sr/⁹⁰Y mixture is acidified to 3 M nitric acid and passed through a series of three Sr• Spec columns in tandem. Essentially all of the strontium present is retained on these columns. Thus, the effluent from the third column contains only ⁹⁰Y, any ⁹⁰Zr daughter formed during the growth period, and any adventitious
20 impurities. This effluent is filtered through a bed of Amberlite XAD-7, then taken to dryness. The residue is taken up in the minimum possible volume of 2 M nitric acid and passed through a column of RE• Spec resin. The

RE•Spec column is defined as a column containing a yttrium-selective extraction chromatographic resin. This resin is prepared with purified Amberchrom CG-71ms resin impregnated with a 1 M solution of octyl(phenyl)-N,N-diiso-butylcarbamoylmethylphosphine oxide (CMPO) in tributylphosphate (TBP). This material, referred to as RE• Spec™ column material (for rare earth specific) is also available commercially from EIChroM Industries. Here, ⁹⁰Y is retained, while any residual ⁹⁰Sr and most impurities are not. Sorbed yttrium is then eluted from the column with a small volume of nitric acid solution. Finally, each of the Sr• Spec columns is stripped of sorbed ⁹⁰ Sr with dilute (~ 0.05 M) nitric acid and the strontium returned to storage to allow for new ⁹⁰Y ingrowth.

It is important to point out that in the processing of the quantities of material typically required for clinical applications, the various stages of the process must be physically isolated from one another to minimize the risk of contamination of the final product with ⁹⁰Sr. In the arrangement now employed at Argonne National Laboratory for the processing of multi-Ci levels of activity, the initial purification of the ⁹⁰Sr is carried out in one containment box, the tandem Sr• Spec column treatment in a second, and the final ⁹⁰Y concentration and purification in a third. Note that clinical applications require that the ⁹⁰Y product be in an aqueous solution containing \leq 0.05 M hydrochloric acid rather than nitric acid, making a cross-over step necessary before final packaging. This is also performed in the third

total strontium loaded onto the first column should not exceed 20% of its capacity.

Experimental Results

Elution behavior of selected elements on a Strontium selective column

5 As already noted, safe clinical use of ^{90}Y requires the reduction of the ^{90}Sr level by a factor of $> 10^6$. In addition, the yttrium must be separated from any trace elements which could interfere in the antibody labeling process. Because this process often involves the chelation of the ^{90}Y with diethylenetriaminepentaacetic acid (DTPA) conjugated to monoclonal
10 antibodies, those elements having chelate formation constants with DTPA within several orders of magnitude of that of yttrium (for which $\log K_f = 22.05$) are of particular concern (See Table 1).

15

20

| Element | Log K_f |
|---------|-----------|
| Fe | 27.5 |
| Ni | 20.0 |
| Cu | 20.5 |
| Zn | 18.0 |
| Y | 22.0 |
| Hg | 27.0 |
| Pb | 18.9 |

Table 2 summarizes the elution behavior of these elements on the Sr• Spec column, along with that of Sr, Y and Zr. Also included are a number of other elements which might be encountered in the fission product solutions used as ⁹⁰Sr feedstock and several elements considered likely to be retained by the Sr•Spec column on the basis of available thermodynamic data for cation interaction with 18-crown-6 based macrocycles. As can be seen, with the exception of barium, mercury, lead and strontium, each of the test elements is essentially unretained by the column. Each, in fact, can be completely eluted with 10 free column volumes (FCV) or less of 3 M HNO₃. Although barium and mercury are more strongly retained by the column, their elution too is virtually complete by the start of strontium elution. It should be noted that the choice of 30 free column volumes of rinsing before Sr stripping is entirely arbitrary. Nearly twice this volume of 3 M HNO₃ [~ 50 FCV] may be passed through the column before strontium breakthrough occurs. Thus, mercury can be completely eliminated from the strontium fraction. Lead is strongly sorbed by the Sr• Spec resin and under identical conditions, distribution ratios of lead between the resin and nitric acid are an order of magnitude or more larger than the corresponding values for strontium. Lead and strontium are readily separated, however, simply by eluting the Sr with ~ 10 FCV of dilute (e.g. 0.05 M) HNO₃ immediately after the 3 M HNO₃ rinses. Under these conditions, lead remains strongly sorbed. This sorbed lead may be removed from the column using dilute (0.1 M) sulfuric acid or ammonium oxalate solution.

TABLE 2
Elution behavior of selected elements of a SR • Spec column*, †

| Element | Fraction eluting (%) by number of free column volumes | | | | | | |
|---------|---|-------|--------|-------|-------|-------|--------|
| | 3M HNO ₃ | | | | | | 0.05 M |
| | 1-5 | 6-10 | 11-15 | 16-20 | 21-25 | 26-30 | |
| Li | 100 | - | - | - | - | - | - |
| Na | 100 | - | - | - | - | - | - |
| Mg | 100 | - | - | - | - | - | - |
| Al | 100 | - | - | - | - | - | - |
| K | 65.8 | 35.3 | - | - | - | - | - |
| Ca | 98.7 | <1.0 | <1.0 | - | - | - | - |
| Mn | 100 | - | - | - | - | - | - |
| Fe | 100 | - | - | - | - | - | - |
| Co | 99.9 | - | - | - | - | - | - |
| Ni | 100 | - | - | - | - | - | - |
| Cu | 100 | 0.2 | - | - | - | - | - |
| Zn | 100 | 0.2 | - | - | - | - | - |
| Rb | 100 | - | - | - | - | - | - |
| Y | 98.5 | 0.1 | - | - | - | - | - |
| Zr | 91.2 | 0.4 | - | - | - | - | - |
| Ag | 15 | 87.6 | - | - | - | - | - |
| Cd | 100 | - | - | - | - | - | - |
| Cs | 100 | - | - | - | - | - | - |
| Ba | <0.4 | <0.4 | 52.9 | 41.6 | 5.9 | 0.7 | - |
| Hg | (6.0) | (5.6) | (18.8) | 45.3 | 19.6 | 14.1 | (6.9) |
| Pb | - | - | - | - | - | - | - |
| Bi | 100 | - | - | - | - | - | - |
| Sr | - | - | - | - | - | - | 99 |

* Sample, 3MHNO₃ containing the indicated elements; V_{sample} = 600 μ L. Column parameters, 1.0 cm³ bed volume; 50-100 μ m particle size resin; flow rate = 1 mL cm⁻² min⁻¹.

† Because of errors inherent in the ICP-AES method used for quantitation, the fractions shown for a given element may not total 100%. Values in parentheses are subject to considerable uncertainty and are intended only as a guide.

That strontium may be efficiently separated from any of the test elements is further demonstrated by the results presented in Table 3, which shows the decontamination factors achievable with the Sr•Spec column. A

5 decontamination factor is the ratio of the concentration of the test element in the feed stock introduced to the column to the concentration of the same element in the fraction of interest, here the strontium fraction. Larger values of this factor indicate a more complete separation of strontium from the test element, this is, more efficient decontamination of the strontium fraction. As can be seen, the decontamination factors observed are typically quite large.

TABLE 3
Decontamination factors for selected elements on a SR• Spec Column*

| Element | Decontamination factor † |
|---------|--------------------------|
| Li | >76 |
| Na | >80,000 ‡ |
| Mg | >700 |
| Al | >100 |
| K | (>75) |
| Ca | >100,000 ‡ |
| Mn | >4700 |
| Fe | >2700 |
| Co | >1100 |
| Ni | >1100 |
| Cu | >3100 |
| Zn | >1300 |
| Rb | >2200 |
| Y | >2600 |
| Zr | >1800 |
| Ru | >40 |
| Rh | >10 |
| Pd | >11 |
| Ag | >2000 |
| Cd | >5600 |
| Cs | >330,000 ‡ |
| Ba | >500 |
| Hg | >(58) |
| Pb | >570 |
| Bi | >470 |

* All conditions as described in Table 2.

† For values bearing a > sign, the measured concentration in the strontium fraction was at or below 3 X the limit of detection for the element (ICP) or with 2σ of the background (radiometric).

‡ Determined radiometrically.

With regard to the preparation of highly pure ^{90}Y , these results are significant for two reasons. First, they indicate that the purity of the ^{90}Sr feed solution (and, therefore, that of the ^{90}Y product) should actually improve with each use of the process. That is, after the initial treatment of a given feed stock solution on the Sr• Spec columns, it is the strip fraction from these columns which becomes the feed for the next run (see Table 2). This fraction is, as these results show, significantly purer than was the original feed. In addition, the strong strontium sorption by the Sr• Spec column indicates that passage of the feed through the column should dramatically reduce the strontium content of the yttrium fraction. Indeed, in experiments using ^{85}Sr tracer, decontamination factors (here, referring to the ratio of ^{85}Sr in the feed to its concentration in the yttrium fraction) of 10^3 or more are typically observed with a single column. Thus, three columns operating in series as shown in Fig. 1 would be expected to yield a total decontamination factor approaching 10^9 .

Elution behavior of selected elements on a RE• Spec column

Table 4 summarizes the elution behavior of each of the test elements on a RE• Spec column. The conditions chosen for this experiment (i.e. 2M HNO_3 as the eluent and < 2-FCV fractions early in the run) were suggested by preliminary studies which showed iron to be the principal impurity encountered throughout the process. As shown in Fig. 2, the apparent selectivity of the RE•Spec resin for Y^{3+} over Fe^{3+} is at a maximum at 0.2-0.5 M HNO_3 . However, at such low acidities, yttrium retention is quite

poor, with < 6 FCV required to reach peak maximum (implying almost immediate breakthrough). Therefore, an acidity was chosen at which yttrium was more strongly retained ($k' > 10$), yet iron retention had not yet become significant. As can be seen from the table, under these conditions, elution of

5 most of the elements, including strontium, is complete in 4 FCV or less. Although iron is somewhat more strongly retained, its elution is nearly 99% complete after only 8 free column volumes. Yttrium breakthrough occurs at 7-8 FCV, with the yttrium concentration in the effluent reaching a maximum at ~ 12 FCV. Its elution is complete by 18 FCV. Several of the test

10 elements, most notably Bi and Zr, are more strongly retained than is yttrium by the RE• Spec column. In each case, however, elution of yttrium is essentially complete before appreciable breakthrough of the element occurs.

TABLE 4
Elution behavior of selected elements of a RE•Spec column*,†

| Element | Portion eluting (%) in each fraction | | | | | | | | | | |
|---------|--------------------------------------|------|------|-----|------|------|------|-----|-----|--------|------|
| | 2M HNO ₃ | | | | | | | | | 0.05 M | |
| | 1 | 2 | 3 | 4 | 5 | 6 | 7 | 8 | 9 | | 10 |
| Li | 79.3 | 27.8 | - | - | - | - | - | - | - | - | - |
| Na | 74.7 | 30.5 | - | - | - | - | - | - | - | - | - |
| Mg | 78.5 | 28.3 | - | - | - | - | - | - | - | - | - |
| Al | 79.4 | 26.8 | - | - | - | - | - | - | - | - | - |
| K | 81.8 | 27.3 | - | - | - | - | - | - | - | - | - |
| Ca | 75.3 | 33.7 | - | - | - | - | - | - | - | - | - |
| Mn | 45.8 | 61.2 | - | - | - | - | - | - | - | - | - |
| Fe | <0.5 | 53.4 | 44.9 | 2.2 | <1.1 | - | - | - | - | - | - |
| Co | 75.0 | 30.3 | - | - | - | - | - | - | - | - | - |
| Ni | 77.3 | 28.2 | - | - | - | - | - | - | - | - | - |
| Cu | 76.4 | 30.1 | - | - | - | - | - | - | - | - | - |
| Zn | 77.2 | 30.2 | - | - | - | - | - | - | - | - | - |
| Rb | 75.9 | 27.2 | - | - | - | - | - | - | - | - | - |
| Sr | 78.8 | 28.9 | - | - | - | - | - | - | - | - | - |
| Zr | 1.4 | 1.1 | 0.3 | - | - | - | - | - | - | - | 76.0 |
| Ag | 82.0 | 17.9 | - | - | - | - | - | - | - | - | - |
| Cd | 72.6 | 34.2 | - | - | - | - | - | - | - | - | - |
| Cs | 74.8 | 27.0 | - | - | - | - | - | - | - | - | - |
| Ba | 79.9 | 27.2 | - | - | - | - | - | - | - | - | - |
| Hg | 47.5 | 51.2 | - | - | - | - | - | - | - | - | - |
| Pb | 63.0 | 41.9 | - | - | - | - | - | - | - | - | - |
| Bi | - | - | - | - | - | - | - | - | - | - | 7.8 |
| Y | - | - | - | 3.8 | 20.8 | 30.3 | 28.9 | 8.3 | 1.8 | - | <1.0 |

* Fractions 1-6 correspond to 1.8 FCV each, fractions 7 and 8 to 3.6 FCV each, and fractions 9 and 10 to 9.0 FCV each.

† All other conditions are as described in Table 2, except that the sample is in 2M HNO₃. Also, ~ 6 x 10⁻³ M oxalic acid is present in the sample solution to solubilize zirconium.

TABLE 5
Decontamination factors for selected elements on a RE• Spec Column*,†

| Element | Decontamination factor † |
|---------|--------------------------|
| Li | >48 |
| Na | >79 |
| Mg | >360 |
| Al | >70 |
| K | (>10) |
| Ca | >600 |
| Mn | >2500 |
| Fe | 180 |
| Co | >770 |
| Ni | >770 |
| Cu | >1700 |
| Zn | >1740 |
| Rb | >580 |
| Sr | >3900 |
| Zr | >1040 |
| Ag | >1390 |
| Cd | >2970 |
| Cs | >1000 |
| Ba | >350 |
| Hg | (>20) |
| Pb | >370 |
| Bi | >220 |

* All conditions as described in Table 4.

† Based on the measured concentrations in fraction (5), the fraction after yttrium breakthrough occurs.

Product purity/radiolabeling efficiency

The radiochemical purity of the ⁹⁰Y product is determined by a combination of the decontamination from ⁹⁰Sr achieved with the series of Sr•

Spec columns and that obtained with the RE• Spec column. The measured overall reduction in the ^{90}Sr level so obtained, a factor $>4.5 \times 10^8$, compares favorably to that observed for previously described generator systems.

5 The chemical purity of the product, that is the absence of appreciable amounts of cations capable of competing with yttrium for binding sites on an antibody, is indicated by the results of antibody labeling studies. An average of 95% incorporation was observed, a value comparable to or higher than results obtained with ^{90}Y produced by other purification schemes.

Radiolytic stability of Sr• Spec and RE• Spec resin

10 Several of the ^{90}Y generator systems previously proposed have been plagued by radiolytic degradation of the resin or extractant used in the separation of ^{90}Y from ^{90}Sr . Earlier work at Argonne has demonstrated that the SREX and TRUEX extraction process solvents, upon which the Sr• Spec and Re•Spec resins used here are based, exhibit reasonably good
15 radiation stability. Nonetheless, because the resins employ modifications of these process solvents and because of the lack of information regarding the radiolytic stability of the support material, the effect of irradiation upon the elution behavior of selected elements on the resins was examined.

20 Table 6 shows the effect of increasing doses of ^{60}Co γ radiation upon the distribution ratios of ^{85}Sr and ^{241}Am (the latter chosen as a convenient gamma-emitting "stand in" for ^{90}Y) between nitric acid solutions and the Sr• Spec and RE Spec resins, respectively. Over the range of absorbed doses considered, americium distribution ratios on RE•Spec are altered only slightly

by irradiation, implying that yttrium elution should not be greatly affected.

This is confirmed by the results presented in Fig. 3, which shows the elution profile for ^{90}Y on RE• Spec columns prepared from both irradiated and

nonirradiated resins. As can be seen, the elution profiles are nearly

5 indistinguishable. This is not entirely unexpected, given the good radiation stability of the corresponding liquid-liquid system, noted above. The stripping of sorbed strontium from Sr• Spec with dilute acid is also largely unaffected by irradiation (Fig. 4). It should be noted, however, that stripping with pure water becomes much more difficult upon radiolysis. As shown in Fig. 5,

10 however, strontium retention is adversely affected at sufficiently high levels of ^{60}Co exposure. Specifically, 3.6×10^7 rads γ irradiation results in a shift in the peak maximum from ~ 97 FCV to 67 FCV. Curiously, the peak width (at half maximum) is considerably narrower with the irradiated resin. As a result, despite the sizeable drop in the elution volume to peak maximum, the
15 breakthrough volume is not significantly lower for the irradiated material.

The reason for this observation remains unclear at present. Additional studies of the radiolytic degradation of Sr• Spec are clearly warranted .

TABLE 6
Effect of radiolysis upon the uptake of strontium and americium by SR•SPEC and RE•Spec resins, respectively*

| Radiation Dose (rads) | D_w † | | | |
|-----------------------|----------|-----------|----------|-----------|
| | SR•SPEC | | RE•Spec | |
| | Low Acid | High Acid | Low Acid | High Acid |
| 0 | 2.34 | 248 | 8.38 | 287 |
| 4.5×10^6 | 2.59 | 254 | 6.32 | 260 |
| 9.0×10^6 | 2.51 | 236 | 6.37 | 265 |
| 1.8×10^7 | 3.32 | 232 | 7.57 | 258 |
| 3.6×10^7 | 3.58 | 249 | 9.47 | 253 |

* Conditions: ~100 mg of indicated resin plus 2.00mL of 0.05M HNO₃ (for all "low acids" measurements), 2 M HNO₃ (RE.Spec, "high acid" measurements), or 3 M HNO₃ (Sr.Spec, "high acid" measurements). ²⁴¹Am or ⁸⁵Sr tracers.
† Calculated from the equation:

$$D_w = \left(\frac{A_0 - A_s}{w} \right) / \left(\frac{A_s}{V_{aq}} \right)$$

where A_0 and A_s represent the activity of the aqueous phase before and after equilibration with the resin, respectively, w is the weight of resin taken, and V_{aq} is the volume of aqueous phase. The number of free column volumes to peak maximum (k') is obtained by dividing D_w by 2.19 and 1.85 for Sr•Spec and RE•Spec, respectively.

⁹⁰Y can be conveniently and efficiently separated from ⁹⁰Sr by extraction chromatography using a combination of Sr-selective and rare earth-selective chromatographic resins. The ⁹⁰Y product obtained is of extremely high chemical and radiochemical purity. No trace elements are present at levels capable of impeding antibody labeling. Strontium-90 decontamination factors approaching 10⁹ are achievable.

Because the ^{90}Sr is removed from the chromatographic columns and stored either as a solution or as a solid between runs and because the extractants used exhibit reasonably good radiolytic stability, rapid deterioration of the resins is not observed, making their frequent replacement unnecessary. Also, because the ^{90}Sr stock is repurified each time that ^{90}Y is recovered, radiolytic and chemical impurities do not accumulate in the stock solution, reducing the frequency with which it must be replaced. Thus, the volume of waste generated by the process is minimized .

The foregoing description of preferred embodiments of the invention has been presented for purposes of illustration and description. It is not intended to be exhaustive or to limit the invention to the precise form disclosed, and obviously many modifications and variations are possible in light of the above teaching. For example, the data from Tables 2 and 4 suggest that other radionuclides of possible interest in nuclear medicine could be prepared using the extraction resins described herein, although modification of the columnar configuration might be necessary. More particularly, the data in Table 2 shows that it is possible to produce ^{82}Rb from ^{82}Sr . The embodiments were chosen and described to best explain the principles of the invention and its practical application and thereby enable others skilled in the art to best utilize the invention in various embodiments and with various modifications as are suited to the particular use contemplated. It is intended that the scope of the invention be defined by the claims appended hereto.

ABSTRACT

An extraction chromatographic method for the preparation of ^{90}Y of high chemical and radiochemical purity is disclosed. After an initial purification of a ^{90}Sr stock solution and a suitable period of ^{90}Y ingrowth, the solution is passed through a series of strontium-selective chromatographic columns, each of which lowers the ^{90}Sr content of the mixture by a factor of about 10^3 . The ^{90}Y remaining is freed from any residual ^{90}Sr , from its ^{90}Zr daughter, and from any remaining impurities by passing the sample through a final column designed to selectively retain yttrium.

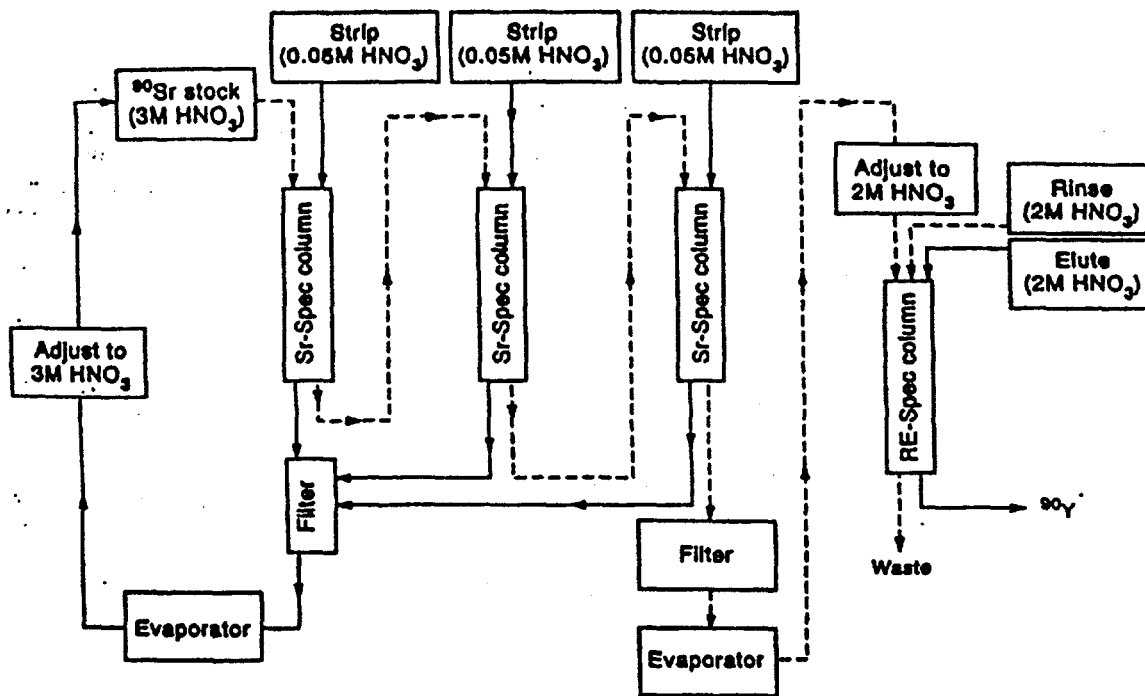


Fig. 1

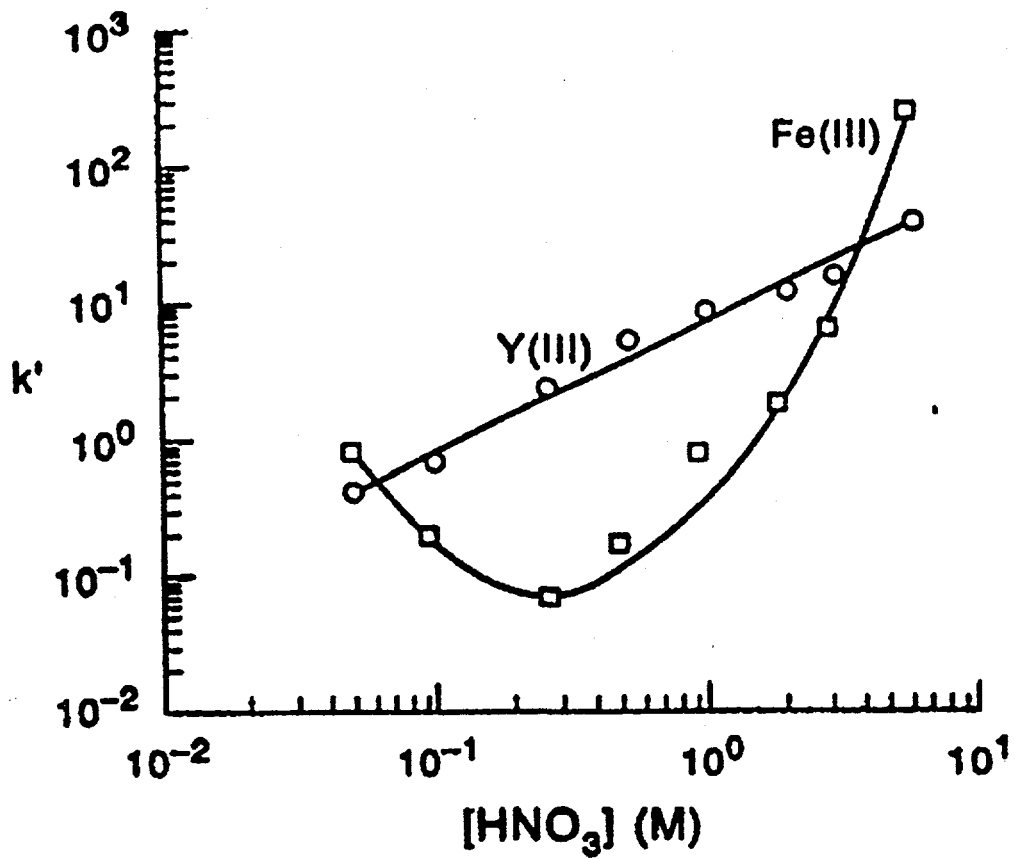


Fig. 2

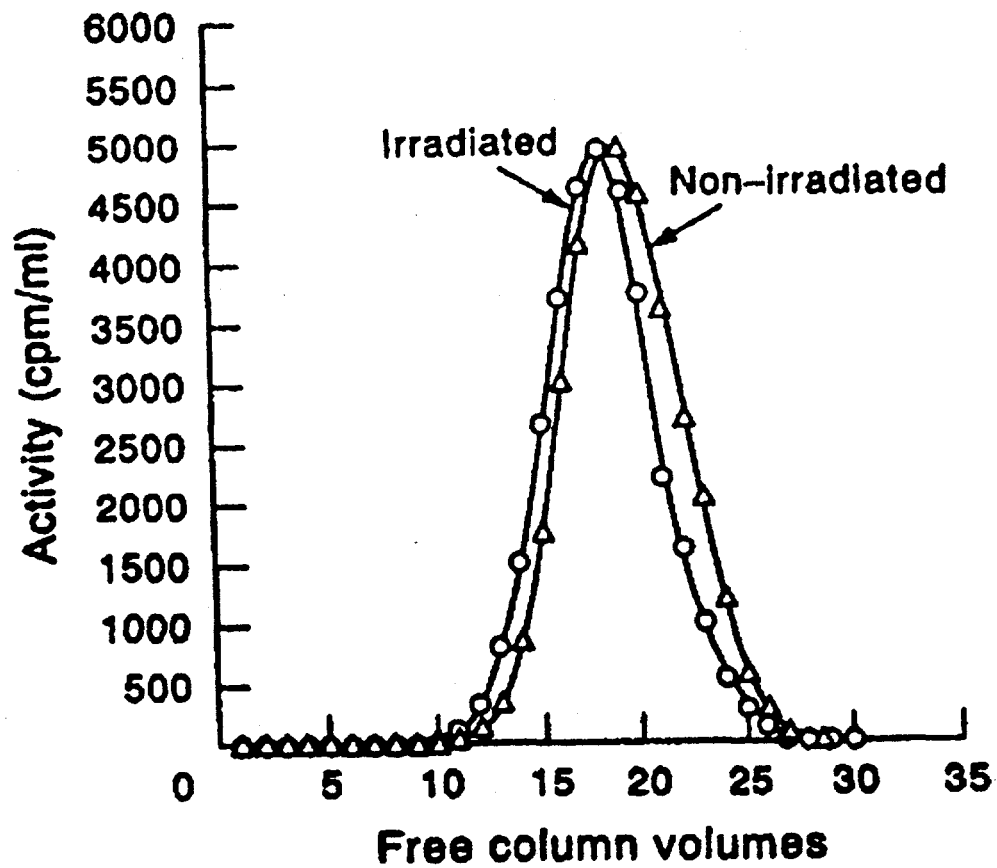


Fig. 3

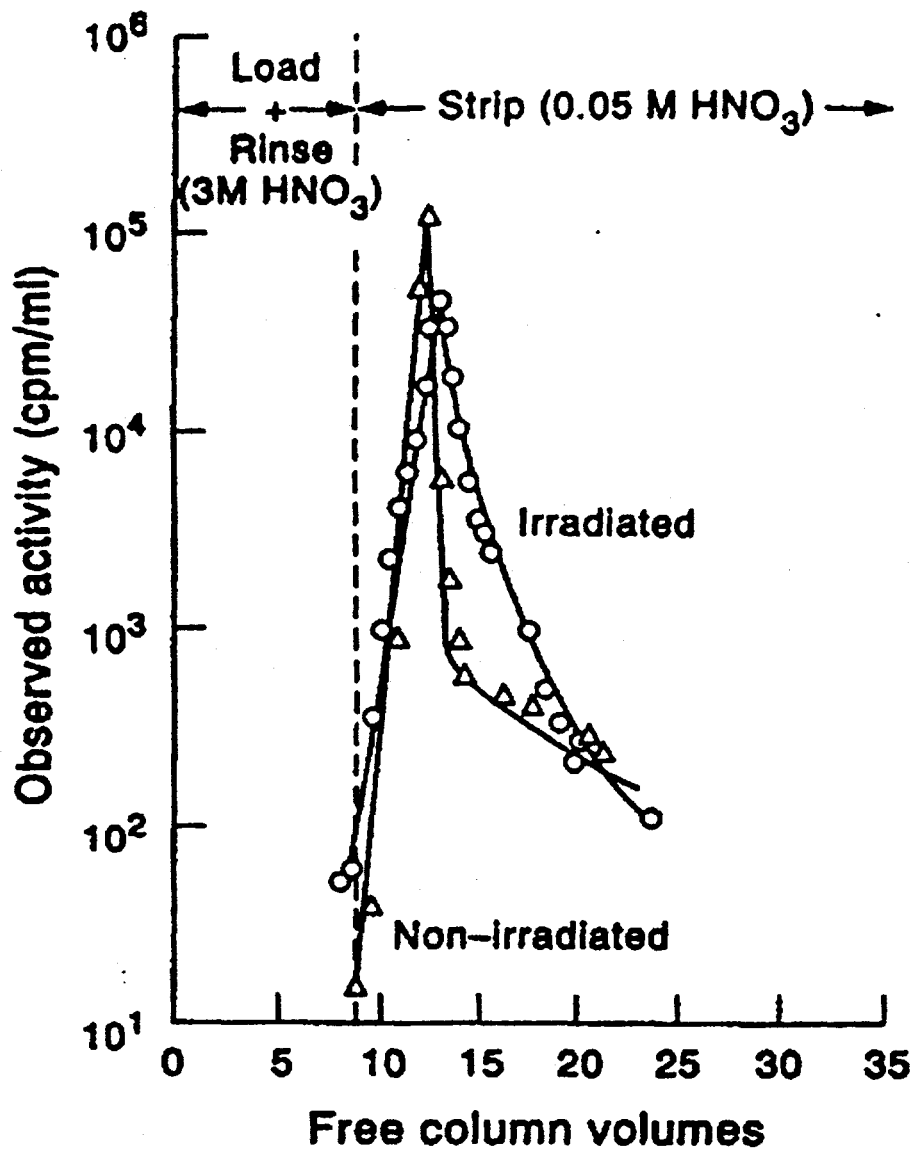


Fig. 4

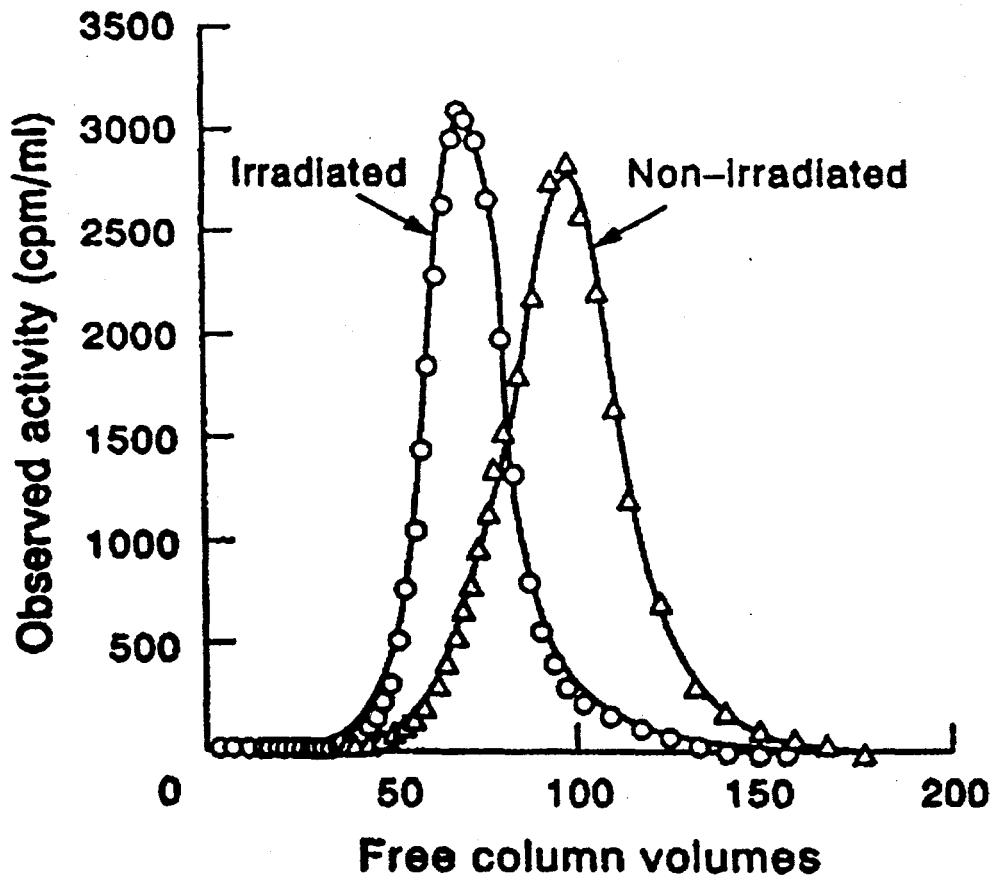


Fig. 5.